



StorAX Nativ

Tubes for stool samples
For in vitro diagnostic use

Intended use

Transport tube without solution for stool transport.

Package contents

25 Stool collection tubes

Storage

The StorAX Nativ tubes are stored at 15°C-30°C. Do not freeze. Do not use the product after the expiration date specified on the package. Once the packaging has been opened, test tubes can be used up to the expiry date if stored properly.

Warnings and Precautions

For professional, single use only. Do not use if evidence of microbial contamination, discoloration, ruptures, or other signs of deterioration; ensure proper product disposal.


Obtaining stool samples

- ① Unfilled StorAX Nativ.
- ② Write the variable data onto the label on the tube.
- ③ Use the spoon to place the freshly obtained stool sample into the transport tube. A spoon-sized sample is sufficient (5 gr. approximately). On suspicion of infestation by intestinal parasites, it is necessary to take at least 3 stool samples on three days in order to isolate eggs and encysted forms.
- ④ Screw on the lid tightly.


1




2

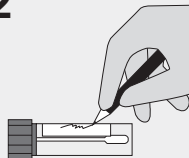

3

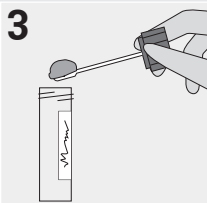

4

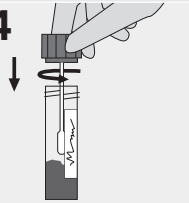

Transporting stool samples with StorAX Nativ

The stool samples can be transported at 15°C-30°C.

Disposal

For any waste produced, each laboratory is responsible for proper disposal in accordance with legal requirements.

Reporting information

A report must be submitted to the responsible authority of the Member State in which the user is located and to the manufacturer in the event of a serious incident in connection with this product.

Abbreviations/Symbols

Symbol	Signification
	Catalogue number
	In vitro diagnostic medical device
	Contains sufficient for <n> tests
	Manufacturer
	Temperature limit
	Use-by date
	Batch code
	Consult instructions for use
	Do not re-use
	CE marking

Change history: Adapted to IVDR specifications

**Axonlab**

Headquarters

Axon Lab AG

Täferstrasse 15
CH-5405 Baden-Dättwil

www.axonlab.com

Affiliates

Austria
Belgium
Croatia
Czech Republic
Germany

info@axonlab.at
info@axonlab.be
info@axonlab.hr
info@axonlab.cz
info@axonlab.de

Luxembourg
Netherlands
Slovenia
Switzerland (de)
Switzerland (fr)

info@axonlab.lu
info@axonlab.nl
info@axonlab.si
info@axonlab.ch
info_f@axonlab.ch

062022 C 03030007